

Tactoset

Injectable Bone Substitute
Small Bone Information Sheet

Tactoset[®] supports cell-mediated regeneration of new bone as material is resorbed^{1,2}

Enhanced flowability & injectability


- Contains hyaluronic acid (HA)
- HA enhances delivery of the calcium phosphate
- May be combined with autogenous bone marrow

Closed mixing system

- Minimizes environmental exposure and enhances ease of use

Augmenting with Tactoset drives surgical confidence


- Increases the density of poor-quality bone caused by cysts or osteoarthritis, resulting in additional suture anchor pull-out strength
- Allows for placement of anchors in the ideal repair location, eliminating the need for extra anchors or complex techniques

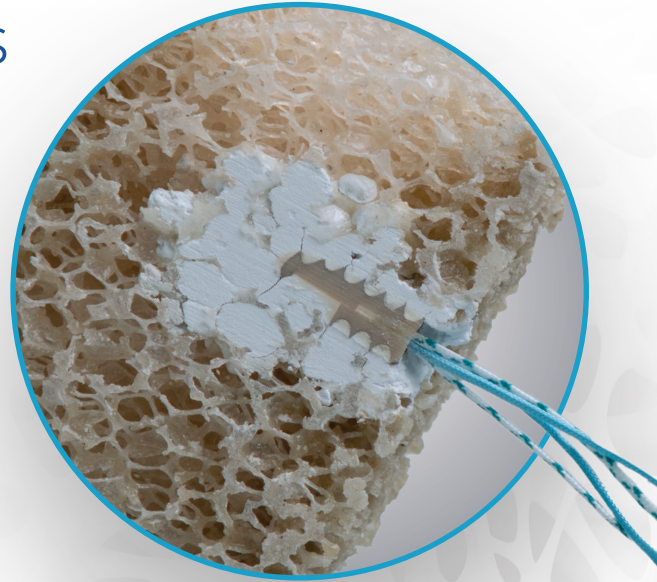

Representation of anchor augmented with Tactoset within trabecular bone

Anchor Pullout Strength


4.5mm Twist PEEK Screw-In Suture Anchor in 10 PCF Solid Foam Block

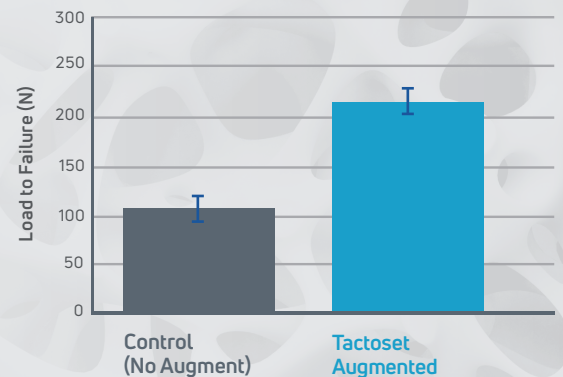

Tactoset has been shown to double soft tissue anchor pullout at body temperature after 24 hours

TACTOSET MATERIAL SET

1. Pre-filled mixing syringe containing powder
2. Glass syringe containing setting solution
3. Four 1 mL Luer lock syringes
4. Female-female Luer lock connector

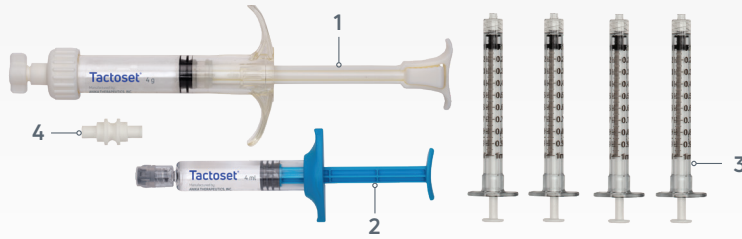

ANIKA SMALL BONE 15GA DELIVERY CANNULAS

1. Drillable outer cannula (with removable drillable stylet)
GAUGE 15 Ga (1.83 mm)
LENGTH 60 mm
2. Push rod



ENHANCED 15GA CANNULA BENEFITS

- Streamlined connector allows multiple pin placements without disruption
- Locks under power
- Sharper cutting angle at the distal tip requires less pressure to get through hard bone
- Minimizes the possibility of overheating

TACTOSET ORDERING INFORMATION

Part Number	Description
1000040	Medium bone 11Ga delivery cannula set including: one drillable outer cannula with removable drillable stylet, one inner cannula for end delivery, one inner cannula for side delivery, one push rod
4000050	Small bone 15Ga delivery cannula set including: one drillable outer cannula with removable drillable stylet, one push rod
6000041	Tactoset material set including: one prefilled mixing syringe containing Tactoset powder, one glass syringe containing Tactoset setting solution, four 1mL Luer lock syringes, one female-female Luer lock connector
6000151	Marrow Cellution Bone Marrow Aspiration Needle set including: one introducer needle with sharp stylet, one aspiration cannula, one 11Ga blunt stylet, one 10 mL syringe

Please note that both the Tactoset material kit item code and one of the delivery cannula item codes are needed for each procedure.

Anika offers a comprehensive portfolio of soft tissue anchors in a broad range of size, material, and suture options. For additional information about Anika's suture anchor portfolio, visit www.anika.com or contact your local sales representative.

For complete product information, please visit: www.anikaifu.com

1. Data on file, Anika Therapeutics, Inc.
2. Tactoset Instructions for Use. Anika Therapeutics, Inc.; AML 500-335.

Anika Therapeutics, Inc.

32 Wiggins Ave., Bedford, MA 01730

1-888-721-1600 | customerservice@anika.com

www.anika.com | Anika. Restore Active Living.™ | Stay Active®

This document and information is intended for markets where regulatory approval has been granted. Anika, Tactoset, Stay Active, and Restore Active Living are trademarks and/or registered trademarks of Anika Therapeutics, Inc. and its affiliates in certain jurisdictions. System designed and manufactured in the USA | Printed in the USA ©2023 Anika Therapeutics, Inc. All rights reserved. AML-900-150 REV. 02