

Integrity™
Implant System
Product Brochure



Strength in Motion.

The Integrity Implant is a synthetic scaffold for rotator cuff repair that provides reliable strength, regenerative biology and a streamlined, precise technique.

Integrity

Implant System

The **Integrity™ Implant** is a porous, knitted scaffold constructed of Anika's **Hyaff®** material, a proven hyaluronic acid technology that supports tissue regeneration and resorbs over time, reinforced with non-absorbable PET (polyethylene terephthalate).

Reliable strength


- Integrity provides higher tensile strength, suture retention, and tear resistance in a thinner, knitted format^{1,2}
- Inherently strong scaffold can be confidently manipulated arthroscopically and enables surgical versatility

Regenerative biology

- Hyaluronic acid scaffold supports regenerative healing through improved cell infiltration, tissue remodeling, and tendon thickening
- Porous, flexible construct knitted using Hyaff, Anika's proven solid esterified hyaluronic acid technology that resorbs over time as tissue remodels

Streamlined technique


- "Lateral first" fixation, unique instrumentation, and simplified surgical technique allow for precise implant placement

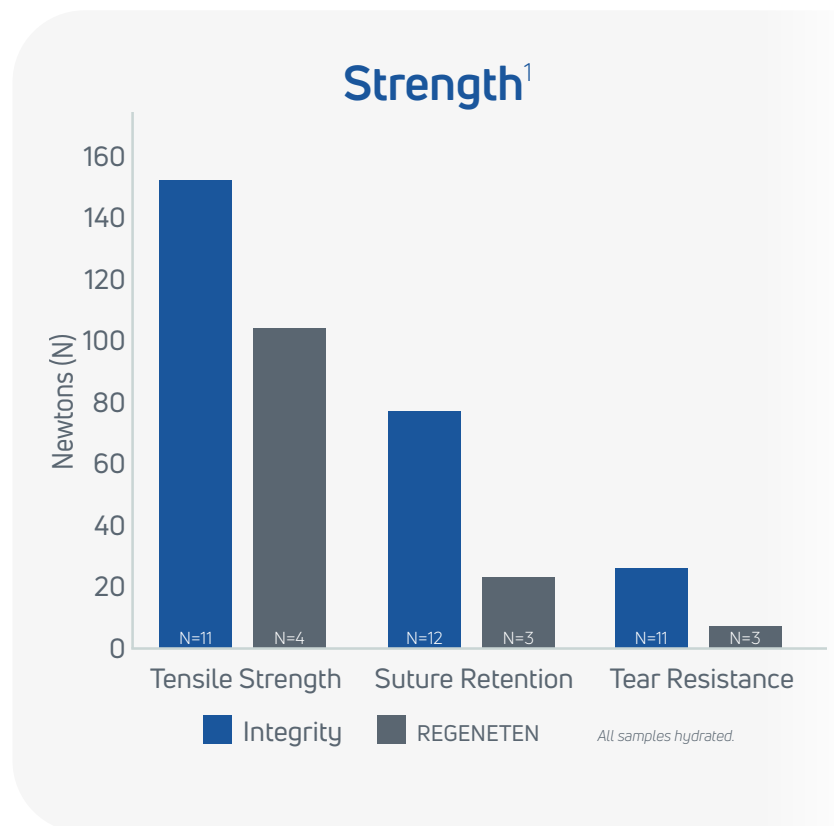

STUDY DATA


A head-to-head animal study was performed comparing the Integrity Implant to REGENETEN.

- As early as 12 weeks post-implantation, fibroblast infiltration and regularly oriented new collagenous tissue formation had occurred within the Integrity repair, demonstrating greater regenerative capacity REGENETEN¹
- At 26 weeks, within the resorbing Integrity structure, new collagenous tissue infiltration forming a new network of tendon tissue had occurred, resulting in nearly 3 times greater thickness in the repaired tendon than REGENETEN¹
- Tendon thickness is thought to improve the local biomechanical environment of the tear by reducing tendon strain, thus optimizing its healing potential³

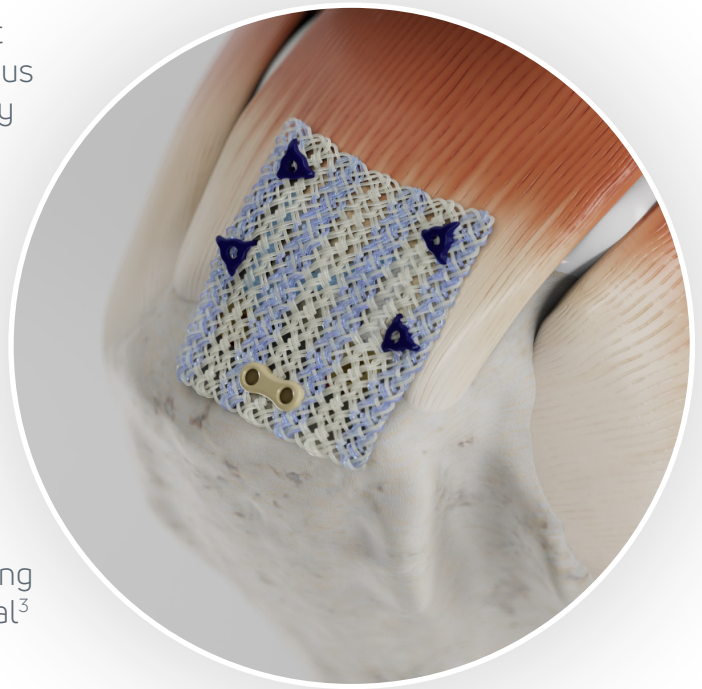


HISTOLOGY


Histology shows Integrity's increased tendon thickness, nearly 3 times greater than REGENETEN.¹

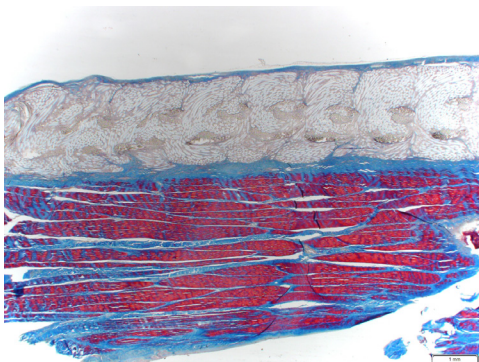

6 weeks

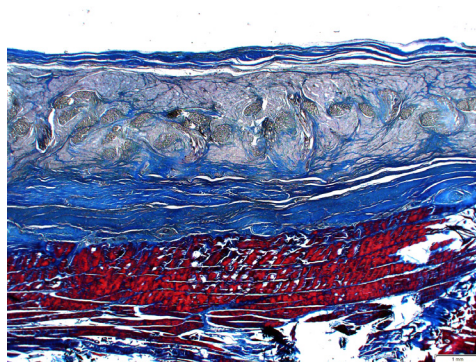

12 weeks

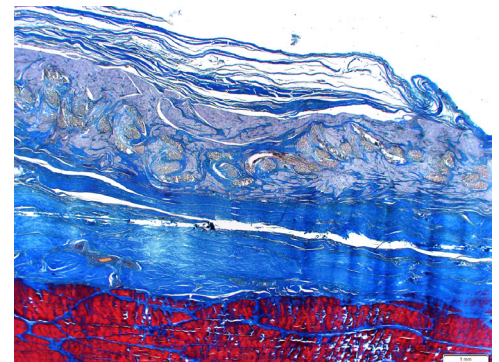

26 weeks

Integrity™ Implant System

INTEGRITY IMPLANT SYSTEM

1. Integrity Implant
2. Implant Fixation Kit, including PEEK Bone Staples and PLGA Tissue Tacks
3. Cannula and Obturator
4. Tissue Tack Delivery Instrument
5. Bone Staple Delivery Instrument
6. Retractable Implant Inserter
7. Deployment Skid



INTEGRITY ORDERING INFORMATION

Part Number	Description
6000100	20x25mm Integrity Implant
6000101	25x30mm Integrity Implant
6000113	25x60mm Integrity Implant
6000114	40x60mm Integrity Implant
6000102	Integrity Implant Fixation Kit, including PEEK Bone Staples and PLGA Tissue Tacks
6000118	Integrity Delivery Instrument, Bone Staple
6000120	Integrity Delivery Instrument, Tissue Tack
6000122	Integrity Cannula/Obturator Kit
6000127	Integrity Retractable Implant Inserter
6000128	Deployment Skid (Reusable Instrument)

1. Data on file, Anika Therapeutics, Inc.
2. Denotes nominal thickness of Integrity Implant
3. Schlegal, T. F., M.D. (2017). Radiologic and clinical evaluation of a bioabsorbable collagen implant to treat partial-thickness tears: A prospective multicenter study. Journal of Shoulder and Elbow Surgery. <https://doi.org/10.1016/j.jse.2017.08.023>

Anika Therapeutics, Inc.
32 Wiggins Ave., Bedford, MA 01730
1-888-721-1600 | customerservice@anika.com
www.anika.com | Anika. Restore Active Living.® | Stay Active.®

This document and information is intended for markets where regulatory approval has been granted.
Anika, Integrity, Hyaff, Stay Active, and Restore Active Living are trademarks of Anika Therapeutics, Inc.
REGENETEN is a trademark of Smith+Nephew.
©2026 Anika Therapeutics, Inc. All rights reserved.
AML-900-557 REV 02

